

# Quick Start App Installation & Activation





## What's in the Box?

### The Perfect Match!

MagHealy comes with an all-new packaging design and what's inside the box is equally impressive. Besides the MagHealy device, you will also find a comprehensive accessories box including an inductive charger, a USB-C charging cable and a manual providing an easy start.

All components are of noble design and built to match each other, allowing MagHealy to always remain front and center, as it deserves.


Wireless **MagHealy** device

**MagHealy** inductive charger

USB charging cable



## Step 1

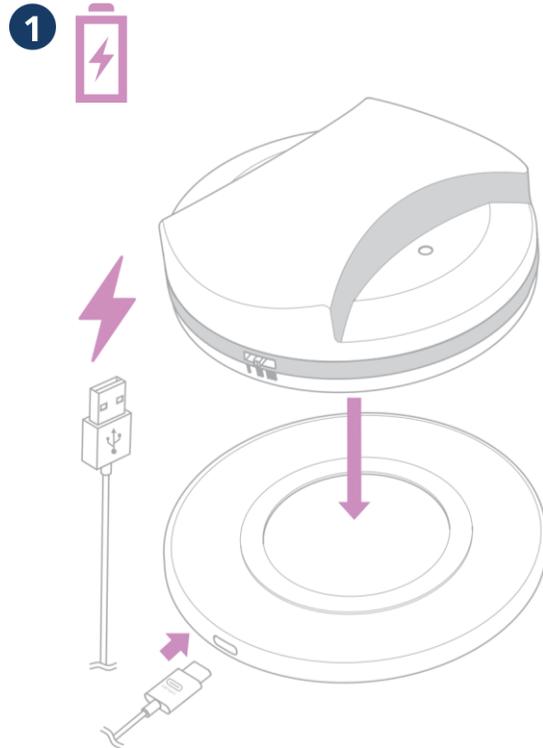

Ensure you have your new MagHealy charged up.

## Step 2



New look




Go to Apple App Store or Google Play Store to search for 'MagHealy' App

- Allows you to configure and control the MagHealy.
- Let you send application programs to the MagHealy.
- Start and stop the application of the MagHealy.
- Connect via Bluetooth to your MagHealy
- You can only connect to one MagHealy at a time.

## Step 3

3 

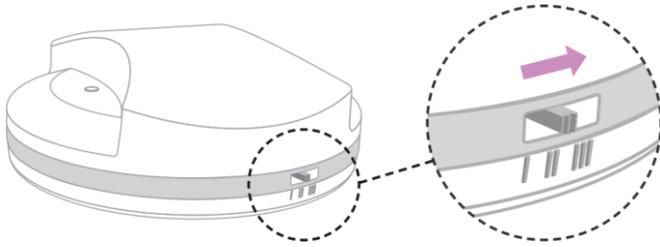


### What do the different positions of the On/Off switch of the MagHealy mean?

In position I, the MagHealy is ready to charge, but the Bluetooth chip and the LED ring are switched off. In position II, the MagHealy is ready to charge, the Bluetooth chip is switched on and you can connect to the Healy App. But the LED ring is switched off. In position III, the MagHealy is ready to charge, and the Bluetooth chip and LED ring are switched on.

Switch on your MagHealy by placing the switch in 2<sup>nd</sup> position (II)

## Step 4

4   
Healy App

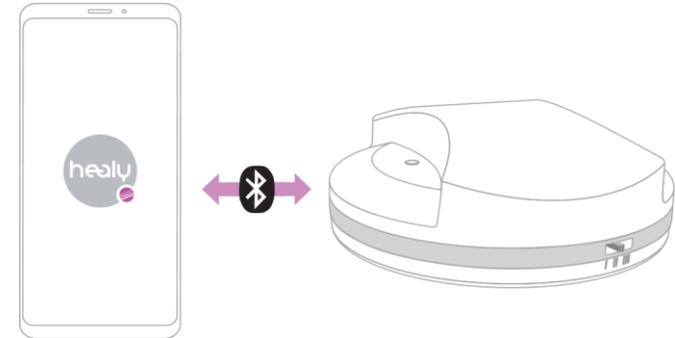

Start the App and connect via Bluetooth to your MagHealy

5




## Step 5: Connecting your device by entering your serial number



Start the App & connect




Scan the QR code & activate

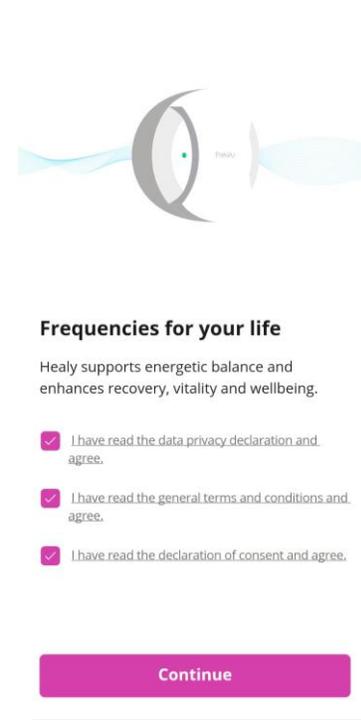

Select the 3 check boxes and click "Continue".



Activate Bluetooth on your phone and click "Search for devices".

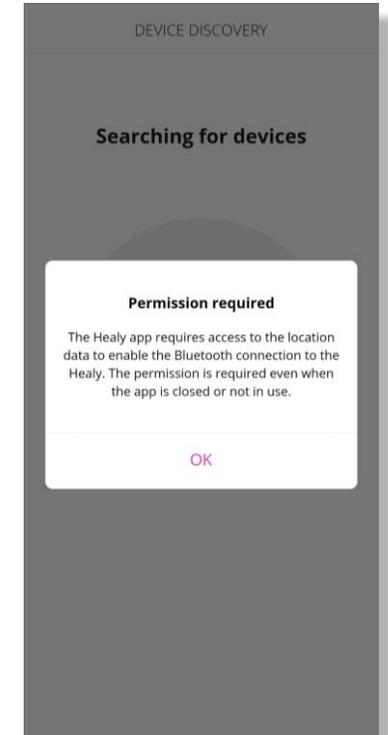

Give permission to location and click "ok".

5



Located on the back of MagHealy:

ESTABLISHING DEVICE CONNECTION


Scan completed, device found!



Scan the code

Enter the code manually

Scan the QR code or enter the serial number manually



Apple permission

CONNECTING TO MAGHEALY



Please enter your device's serial number

03 : 01 : 01 : 01 : 01 : 9D

Register

Click "Register" after entering the serial number


CONNECTING TO MAGHEALY

Please wait



Please wait for your MagHealy to connect.

## Step 6

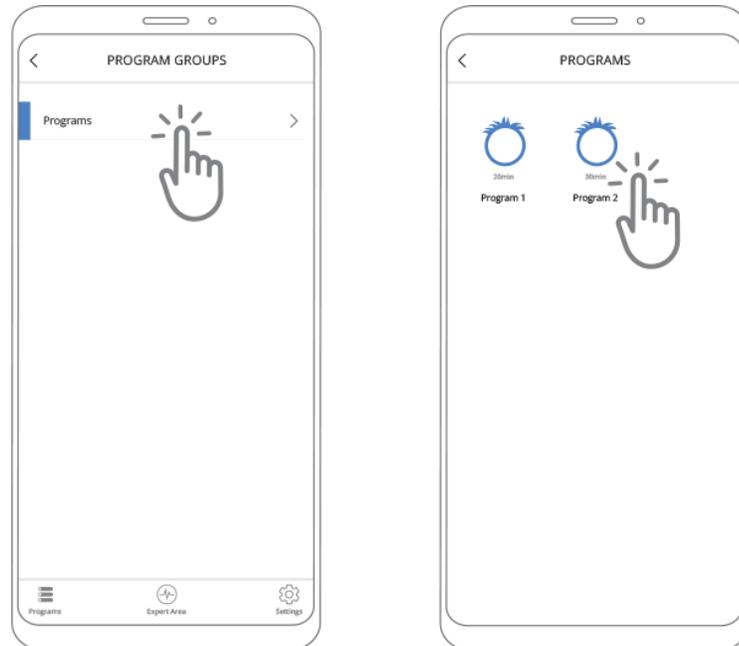6 

Once it's set up, you can now  
start to run your programs

